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Translation of the original instructions
Page 56 – 103

Certifications

Contents

Introduction

About this document.....	56
Typographic conventions.....	56
Copyright protection.....	56
List of abbreviations.....	57

General information

Notes on the product	58
Directives applied, product certification.....	58
California Residents.....	58
Copyright notice	58
Residual risk.....	58
Intended use.....	58
Reasonably foreseeable misuse	59
Transportation.....	59
Storage.....	59
Acclimatization	59
Permissible ambient conditions	59

Safety

General safety instructions.....	60
Electrical safety.....	60
Operational safety	60
Occupational safety	60
Personal protective equipment (PPE).....	61
Environmental protection.....	61
Biohazard.....	61
Special hygiene measures for the use of laboratory equipment in food, cosmetics and pharmaceutical production	61
General measures.....	61
Device-specific measures	62
Other regulations.....	62

Device description

Design.....	63
Operating elements	64
Key combinations	64
Functionality.....	65

Display and user interface.....	65
Display icons	66
Signal sounds	68

Commissioning

Installation and connection	69
Table top version	69
Attachment to Hei-VAP Rotary Evaporator	69
Electrical connection	71
Prepare wall power supply plug	71
Connect to mains	71
Vacuum connection.....	72

Operation

First steps	74
Handling	75
Switch-on/-off controller.....	76
Select language and pressure unit.....	77
Select operation mode	78
Pump down	78
Vac control	78
Program.....	78
Start controlling	79
Venting	79
Change operation mode.....	80
Display graphic (curve)	81
Set vacuum	83
Stopping the controller	83
Operation menus.....	84
Pump down	84
Vac control	85
Program.....	86
Configuration	87
Default settings	87
Autostart.....	88
Display	88
Sensors.....	88
RS-232	89
Function.....	89
Program functions	89

Troubleshooting

Troubleshooting.....	90
Error display	90
Error messages	91
Error messages external components	94
Reset	95
Auto reset	95
Active reset	95

Appendix

Technical data	96
Wetted parts.....	97
Interface commands	97
Pin assignment RS232 interface	97
Scope of delivery.....	98
Accessories	98
Device service.....	99
General cleaning instructions.....	99
Clean venting valve.....	99
Clean internal sensor.....	99
Sensor adjustment	100
Adjustment at atmospheric pressure	100
Adjustment under vacuum	100
Adjustment at a reference pressure	101
Repairs	101
Maintenance.....	101
Disposal.....	102
Contact information Heidolph international.....	102
Warranty Statement	102
Certificate of decontamination.....	103

About this document

These operating instructions describe all the functions and the operation of the Hei-VAC Control vacuum controllers.

The operating instructions are an integral part of the described device!

Typographic conventions

Standardized symbols, signal words and highlighting elements are used in this document to warn of hazards and to identify important information or rather special text contents.

Symbol	Signal word / explanation
	<p>Warning symbols in combination with a signal word indicate dangers:</p> <p>DANGER Indicates an immediate dangerous situation. Failure to observe the indications will result in death or serious injury.</p> <p>WARNING Indicates a potential danger. Failure to observe the indications will result in serious injuries.</p> <p>CAUTION Indicates a potential hazard which, if not avoided, may result in damage to property and minor to moderate injuries.</p>
	<p>Mandatory signs are used to indicate important and useful information on handling a product.</p> <p>This information is used to ensure operational safety and to maintain the value of the product.</p>
[GUI]	<p>Parameter designations, display texts, and device labels are typographically highlighted in the text and in tables to facilitate their assignment on the device.</p>
→	<p>The arrow indicates specific instructions to be followed to ensure operational safety when handling the product.</p>

Copyright protection

This document is protected by copyright and is intended for use by the purchaser of the product only.

No transfer to third parties, reproduction in any form, including excerpts, and by any means, as well as utilization and/or disclosure of the contents is permitted without the prior written consent of Heidolph Instruments GmbH & Co. KG. Any violation is subject to compensation for damage.

List of abbreviations

abs.	absolute
ATM	Standard Atmospheric Pressure
Hei-VAC Control	Vakuüm controller, Controller
di (di)	Interior diameter
DN	Nominal diameter
EK	Emission condenser
EX*	Outlet
FPM	Fluorelastomer
Gr.	Size
hh:mm:ss	Time settings in hour/minute/second
hPa	Pressure unit, hectopascal (1 hPa = 1 mbar = 0.75 Torr)
IN*	Inlet
KF	Small flange
Max	Maximum value
mbar	Pressure unit, millibar (1 mbar = 1 hPa = 0.75 Torr)
Min	Minimum value
min	Minute
PA	Polyamide
PBT	Polybutylene terephthalate
PE	Polyethylene
respon.	responsible, supervising Specialist
RAN-No	Return Merchandise Authorization number
SW	Wrench size (tool)
Torr	Pressure unit (1 Torr = 1.33 mbar = 1.33 hPa)
VAC	Vacuum
VMS	Vakuüm-Management-System

Notes on the product

Directives applied, product certification

CE marking

The device meets all requirements of the following directives:

- Low Voltage Directive 2014/35/EU
- Electromagnetic Compatibility Directive 2014/30/EU

California Residents

Important information for California residents regarding Prop 65. Please visit www.P65Warnings.ca.gov for more information.

Copyright notice

The software implemented in this product is protected by copyright laws. The rights holder is Heidolph Instruments GmbH & Co. KG, Walpersdorfer Straße 12, 91126 Schwabach/Germany. Any open source software components in this software are excluded from our copyright. Further information is available in the service area on our website www.heidolph-instruments.com.

Violations of copyright (e.g. unauthorized use or modification of the software) may result in civil claims (e.g. omission, compensation for damages) and/or entail criminal penalties within the meaning of the legislation of the destination country.

Residual risk

The device has been designed and manufactured in accordance with the state-of-the-art standards at the time of development and the recognized safety regulations. During mounting and use, as well as during maintenance, repair and cleaning work, there are nevertheless certain residual risks associated with the described device.

These are identified and described at the appropriate points in this document.

Intended use

The Hei-VAC Control is a laboratory instrument, used to measure and/or control vacuum in therefore intended plants.

Any other use of this device is not considered as intended!

Due to its design, the device in its delivery condition may only be used in analytical processes or in laboratory-like conditions in the food, cosmetics, and pharmaceutical industries as well as other comparable industries that manufacture products intended for consumption by humans or animals, or for use on humans or animals.

All risks resulting from improper use are borne solely by the operator.

The device may only be operated by authorized and instructed personnel. Training and qualification of the operating personnel as well as ensuring that the device is handled responsibly are the sole responsibility of the operator!

Reasonably foreseeable misuse

Additional measures may be necessary, and/or specific directives and safety regulations may have to be observed for use under conditions or for purposes deviating from the intended use (see, for example, section „Special hygiene measures for the use of laboratory equipment in food, cosmetics and pharmaceutical production“ on page 61). Corresponding requirements must be evaluated and implemented by the operator in each individual case.

Compliance with and implementation of all relevant directives and safety measures for the respective field of application is the sole responsibility of the operator.

Transportation

Keep the original packaging in a dry and protected place for later use.

Storage

Always store the device in its original packaging. To protect against damage and excessive material aging, store the device in an environment that is as dry, temperature-stable and dust-free as possible.

Recommended storage conditions:

- -10 °C – 60 °C at 30 % – 85 % rel. humidity, avoid condensation

Acclimatization

After each transportation and after storage under critical climatic conditions (e.g. high temperature difference between inside and outside), allow the device to acclimatize at room temperature for at least two hours to prevent possible damage due to condensation before commissioning it in the place of use.

If necessary, extend the acclimatization phase if the temperature differences are very high.

Permissible ambient conditions

The device is designed for indoor use only. Permissible ambient conditions for operation:

- 10 °C – 40 °C at 30 % – 85 % rel. humidity, avoid condensation
- Installation altitude: 3,000 m asl

When used in corrosive atmospheres, the service life of the device may be reduced depending on the concentration, duration and frequency of exposure.

The device is **NOT** suitable for outdoor use!
The device is **NOT** suitable for use in potentially explosive areas!

General safety instructions

- Before commissioning and using the device, familiarize yourself with all the safety regulations and occupational safety guidelines applicable at the place of use and observe them at all times.
- Only operate the device if it is in faultless technical condition. In particular, ensure that there is no visible damage on the device itself, the power supply cord and, where applicable, on connected devices or the supply connections.
- If there is missing or misleading information on the device or regarding occupational safety, contact the responsible safety specialist or our technical service.
- Only use the device in accordance with the regulations on intended use („Intended use“ on page 58).

Electrical safety

- Before connecting the device to the power supply, ensure that the voltage indicated on the rating plate matches the specifications of the local power utility company.
- Ensure that the power supply circuit provided is protected by means of a residual-current device (RCD).
- Only use the short-circuit-proof multi-voltage power supply supplied with the device.
- The device must only be supplied with power from a properly grounded mains socket-outlet.
- Have repairs and/or maintenance work on the device carried out exclusively by an authorized electrician or by the technical service department of Heidolph Instruments.
- Always disconnect the device from the power supply system before carrying out any maintenance, cleaning or repair work.

Operational safety

- Do not make any unauthorized changes or modifications to the device!
- Only use genuine spare parts and accessories, or those expressly approved by the manufacturer!
- Rectify malfunctions or faults on the device immediately.
- Switch off and disconnect the device from the power supply, preventing reconnection, if it is not possible to eliminate the malfunction or rectify the fault immediately.
- Observe all other applicable regulations such as laboratory and workplace guidelines, recognized safety technology rules and special local regulations.

Occupational safety

- Always use the prescribed personal protective equipment (PPE) such as protective clothing, safety goggles, protective gloves, safety shoes, etc.
- Depending on the running process, explosive mixtures can form in plants or other hazardous situations could result. Never operate critical processes unattended!

Personal protective equipment (PPE)

The operating company must determine and provide the necessary PPE, depending on the respective application and the media and chemicals used.

The corresponding instruction of the personnel is solely within the operating company's responsibility.

Environmental protection

When processing environmentally hazardous substances, take appropriate measures to avoid hazards to the environment.

The evaluation of corresponding measures such as the marking of a hazardous area, their implementation, and the training of the relevant personnel is the sole responsibility of the operator!

Biohazard

When processing biohazardous substances, take appropriate measures to prevent hazards to persons and the environment, including:

- Instruction of the personnel regarding the necessary safety measures.
- Provision of personal protective equipment (PPE) and instruction of the personnel in its use.
- Marking the device with a biohazard warning symbol.

The evaluation of corresponding measures such as the marking of a hazardous area, their implementation, and the training of the relevant personnel is the sole responsibility of the operator!

Special hygiene measures for the use of laboratory equipment in food, cosmetics and pharmaceutical production

When laboratory equipment is used in the production processes of the food, cosmetics or pharmaceutical industry, special hygiene measures must be taken by the user to avoid sample contamination and to minimize any risk to humans and the environment as far as possible.

General measures

- Ensure a clean working and storage environment when handling substances and materials.
- Train all employees in the field of occupational hygiene, document all training measures and check the implementation of all required hygiene measures regularly during operation.
- Use a hygiene control concept such as HACCP (Hazard Analysis and Critical Control Points). The HACCP comprises the following criteria:
 - Hazard analysis
 - Identification of critical control points
 - Definition of critical limit values
 - Establishment of a system for monitoring and controlling critical control points (CCP)
 - Corrective actions for uncontrollable CCP
 - Establishment of a system to verify the implementation of all HACCP measures
 - Establishment of a system for documenting all associated procedures and records

Device-specific measures

- Regularly clean components that come into contact with the product in the autoclave (if available or possible) or chemically (e.g. with ethanol) to sterilize all surfaces.
- Make sure that even products that are intended for single use only are of sufficient purity.
- Do not use open containers.
- Avoid contamination caused by careless handling of contaminated vessels, apparatus or aids.

Contact information

For further information, please contact our after sales service at any time.

Phone: +49-9122-9920-0

Email: sales@heidolph.de

Other regulations

In addition to the notes and instructions in this document, observe all other applicable regulations such as laboratory and workplace guidelines, hazardous substances ordinances, recognized rules of safety engineering and occupational medicine as well as particular local regulations!

Noncompliance will invalidate any warranty claims against Heidolph Instruments.

The operator is solely liable for all damage resulting from unauthorized changes or modifications to the device, from the use of unauthorized or non-genuine spare parts and accessories, or from disregarding the safety instructions and hazard warnings or the manufacturer's instructions!

Design

Legend

1	Chemically resistant plastic housing
2	LC display
3	Control panel
4	Selection knob
5	Ports for Heidolph vacuum valve
6	
7	Mains connection
8	Serial port RS 232 C (Sub-D)
9	Hose nozzle, vacuum connection
10	Venting tube, connection for external venting, e. g., inert gas
11	Rating plate
12	Spring clip as fixation for built-in version or plastic foot for table top version

Operating elements

Press [Selection knob]

- Call up menu or function.
- Exit menu or function.
- Confirm and store setting or selection.
- Step to next menu, function or content.
- Call up Configuration menu while booting.

Turn [Selection knob]

- Select menu, navigation function.
- Value setting – increase/decrease.
- Adapt set vacuum (in mode Vac control).

[On/Off]

- Switch on/off controller.

[Start/Stop]

- Start/stop vacuum control.
- Confirm completed program when clock icon blinks.
- Confirm error and status indications.

[System venting]

- Keystroke < 2 s = momentarily venting, control continues.
- Keystroke > 2 s = venting to atmospheric pressure (max. 1050 mbar/787 Torr.), control stops.
- Keystroke while venting = venting stops.

[Mode] – Select operation mode

- With stopped operation: Mode menu for selecting the operation mode.

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- During running operation: To switch from Pump down to Vac control and further to Auto mode.
- During running operation: To switch between Auto mode and Vac control.

Key combinations

Menus and functions that are not intended for everyday use, can only be accessed through button combinations (operating elements see table above):

- Press and hold [Selection knob], press [On/off] (when controller is switched off)
 - Call up [Language selection] menu
 - Call up [Pressure unit] menu
- Press and hold [VENT], press [On/off] (when controller is switched off)
 - Call up [Function] menu
- Press and hold [Mode], press [Selection knob] ([Function] menu only)
 - Enable Vacuum address assignment (frame marking).
 - Confirm parameter transfer of selection Vario init.
- Press, hold and turn [Selection knob]
 - Quick adaption set vacuum (in mode Vac control)
 - Quick adaption motor speed high/low in [Pump down] mode – only in combination with VARIO®

Functionality

The controller is designed for applications requiring controlled vacuum. The controller has a two-point control mode to switch an in-line isolation valve.

The controller is freely programmable. Up to 10 programs can be stored in the controller memory. Each program also offers up to 10 program steps (time/pressure) plus control functions, such as: venting, pump down and ramp function.

The controller enables the measurement of relative pressure with regard to a reference sensor (VSK 3000). The controller manages vacuum processes by controlling vacuum pumps, in-line isolation- and/or air admittance valves. It controls process vacuum, cooling water and venting to demand.


Valves and/or vacuum pumps are necessary to operate the controller. Without those components the controller can only be used as vacuum measurement device.

- Ceramic diaphragm vacuum sensor and venting valve are already integrated into the controller.
- The ceramic vacuum sensor is chemically highly resistant, measures accurately and is gas-type independent.
- While booting the controller checks for current configuration of connected components.
- Operating elements are the selection knob, buttons on the control panel and full text menus on the display.

Display and user interface


After booting the pressure display appears, including Bar graphic and preset operation mode:

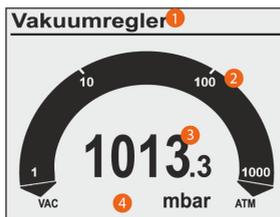Legend

1	Status bar (title bar) Mode: Pump down, Vac control, Program Process time: hh:mm:ss (only displayed with running process)
2	Bar graph: Graphical display of actual pressure
3	NUmeric value: Actual pressure = digital pressure display
4	VAC: Vacuum mbar: Pressure unit according to pre-setting (mbar, Torr, hPa) ATM: Atmospheric pressure

Display Icons

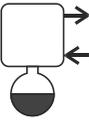Symbol	Meaning
	Vacuum control is running (animation)
00:00:00	Process time; runtime vacuum control (hh:mm:ss)
	Pump is running; in combination with percentage sign = motor speed (only for VARIO systems)
	Venting* valve is active, i. e. open (VENT); Flashing cycle: continous venting switched on.
	Coolant valve switched on, open
	In-line valve switched on, open
	Emission condenser (Peltronic) connected
	Level sensor activated (only when level sensor is connected)
	Pump down – continous pumping
	Pump down: lower pressure limit reached VACUU-LAN: pump down to set pressure Vac control: for 2-point control – pump down to set pressure

Device description

VACUU-LAN: pressure increase to switch on pressure
Vac control: preset maximum exceeded

VARIO control: pump down to set point.
Auto mode: pump down and boiling point detection within the preset time interval regarding changing process conditions.

VARIO control: reaches and tracks boiling point. The next program step starts when the programmed pressure has been reached or the preset time has elapsed.

2-point control: pressure in hysteresis, pump switched on

2-point control: pressure in hysteresis, pump switched off

Turbo mode switched on (for VARIO® pump in combination with turbo-molecular pump)

Clock
- Program completed*
- Mode VACUU-LAN: delay time elapses

Lock – operation locked

HI mode for Pump down = optimum speed for the respective pressure.

Percentage value for Pump down motor speed.

Set value for Vac control.

Flashing: warning!

Signal sounds

Setting Sound On in menu Configuration/Display is required to hear the audio signals..

1x 

Short beep for each keystroke.

2x 

Audio warning for error indication.

In short intervals a number of warning beeps are to be heard. This Audio warning is active until error clearance or reset.

Error messages are indicated by differing numbers of beeps (audio warning).

For the list of possible warning beeps see chapter Fault indication

Installation and connection

CAUTION

- Measuring error due to an obstructed vacuum line. Prevent overpressure > 1060 mbar (> 795 Torr) in the piping system
- Condensate can falsify the measurement. Position the vacuum hose in such a way that condensate cannot flow towards the controller and its vacuum sensor. No liquid should accumulate inside the vacuum hose. Install vacuum hoses in such a way that condensate cannot flow into the controller.
- Particles, liquids or dust may not enter the controller. Install a separator or filter at the intake of the system. Appropriate filters are for example chemically resistant, and resistant to clogging.
- Installation and operation in areas where explosive atmospheres can occur is not allowed.
- The use of gas ballast or the operation of venting valves is only permitted if thereby explosive atmospheres normally do not occur in the interior of the equipment or, if they do occur, are likely to do so only infrequently and for a short period. If necessary vent with inert gas.

Table top version

The table top or bench-top type controller can be installed and connected directly on top of the work bench or on laboratory table.

The table top version is supplied with a hose nozzle. The hose nozzle should be positioned in a way that the connected vacuum hose cannot kink.

Attachment to Hei-VAP Rotary Evaporator

The Hei-VAC Control can be attached to Hei-VAP Rotary Evaporator directly.

- Screw the pre-assembled metal plate that carries the vacuum valve to the back of the evaporator with two screws:

Commissioning

- Attach the support rod to the evaporator and attach the controller to the top end of the support rod, so that the support rod locks into the gap on the back of the controller.
- Connect the power supply to the controller [1].
- Connect the cable of the vacuum valve to the controller [2].

- Attach the hoses:
 - From the upper connector of the vacuum valve to the top of the condensate cooler [3].
 - From the lower connector of the vacuum valve to the vacuum source [4].
 - From the connector for measuring and venting of the controller to the second connector at the top of the condensate cooler [5].

Electrical connection

Wall power supply kit: short-circuit-proof multi-voltage power supply with integrated overload protection and changeable mains plugs:

Prepare wall power supply plug

- Take the wall power supply kit out of the packaging.
- Select the mains plug that fits to your mains socket.
- Connect the mains plug to the metal contacts of the wall power supply plug.
- Slide the mains plug until it locks.

Remove mains plug

- Press the locking knob on top of the wall power supply plug.
- Remove the mains plug.
 - Another mains plug can be fixed.

CAUTION

Please install the power supply line in such a way, that no damage be caused to the cable due to sharp edges, chemicals or hot surfaces.

Connect to mains

- Plug female connection of the power supply cable into mains connection of the controller.
- Plug the wall power supply into the mains socket.
 - Green LED at wall power supply plug glows.

Vacuum connection

CAUTION

Flexible vacuum hoses can contract because of evacuation.

- Fix vacuum hose at the connections.
- Fix connected components.
- Measure and trim the vacuum hose to a length that cares for the maximum shrinkage.

Possible damages to parts which are in contact with process media.

Residuals of aggressive or condensing media can cause damages to the controller or its inner parts.

- Prevent that damaging process media can get into the controller.

Filters will compromise measurement and control.


Connect vacuum line

- Connect the vacuum line gas-tight to the vacuum port of the controller.
- Only use a vacuum hose that is sufficient for the purpose and which provides enough stability.
- Use hose tubes as short as possible.
- Maximum admissible pressure at vacuum sensor: 1,5 bar/ 750 Torr (absolute).
- Observe the maximum measuring limit of the controller, approximately 1060 mbar (795 Torr).

Connection examples


Depending on design and installation the controller provides several options for connection to the vacuum system.

- Flexible caoutchouc hose, directly plugged on the hose nozzle:

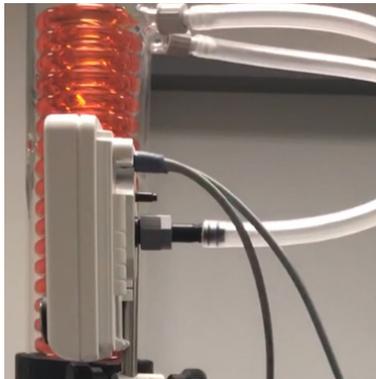Commissioning


→ Vacuum hose made of PTFE, plugged on hose nipple, fixed with union nut:

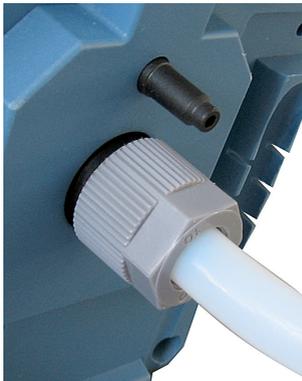With built-in controller the vacuum port is not visible. The vacuum hose may not be kinked.

Make sure there is sufficient space inside the housing, or use a stable, curved hose nozzle for connection.

First steps

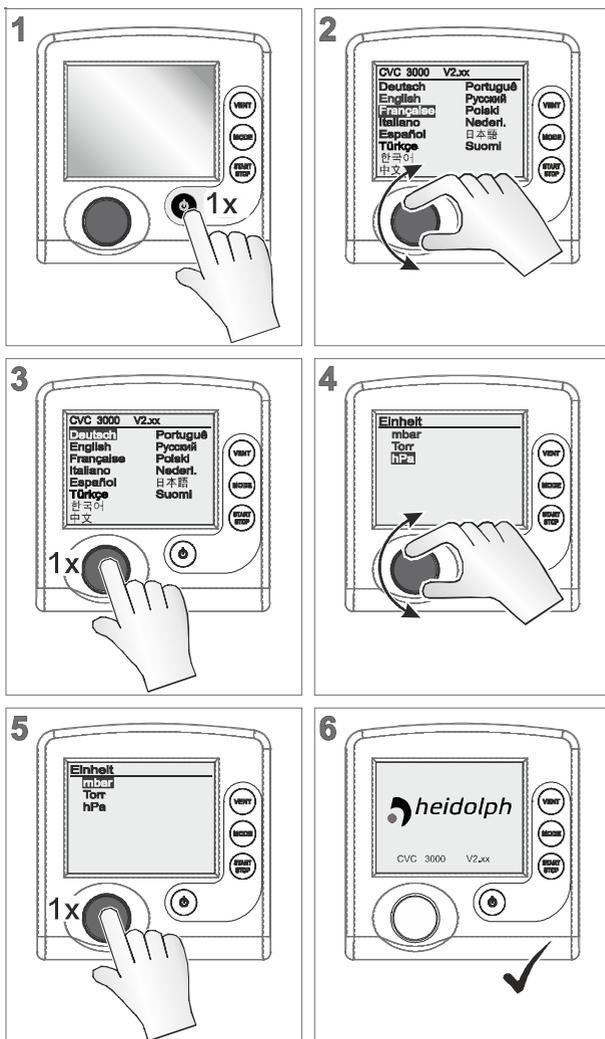
Handling

The handling of the controller is menu-driven. Menus are accessed via push buttons on the control panel or via button combinations. Use the selection knob to select function or menu.

Navigation

- Turn selection knob to select a menu by shifting the bar marking: bar marking up/down:

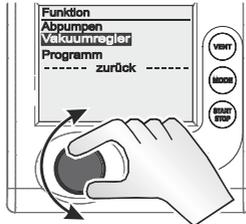
Sub-menus

- Sub-menus are highlighted with points:


Selection

- Press selection knob to confirm selection:

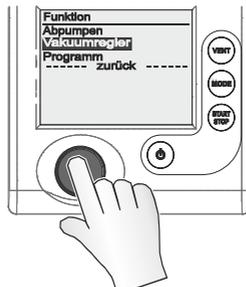Input (data entry)

- Changeable values are positioned on the right side in the display:

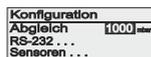- Text on the right side accords to content selection like in a drop down list (Exception: menu Program, in this menu data and value are editable):


Example: enable entry and edit

- Select the required line and press selection knob.
 - Marking jumps to the input field.
 - Input/Content selection enabled.
- Turn the selection knob.
 - Value/Content changes.
- Adapt the numeric value within the specified min/max range or select the required function out of the available content.
- Confirm input/selection by pressing the selection knob.
 - Value is stored or
 - selected function starts


Back (return)

- Place the bar marking on line back and press the selection knob to return to previous menu, display or to pressure display:

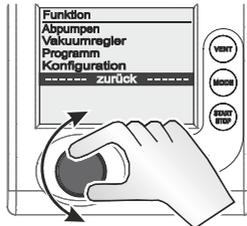
In sub-menu Sensors the display returns to previous menu only after the selection of a sensor.

Switch-on/-off controller

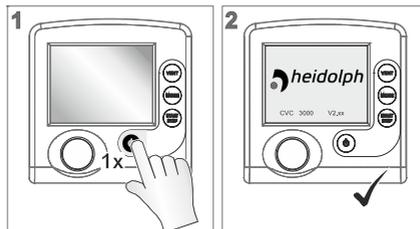
- Switch on by pressing on/off
 - Initial screen: company logo and firmware version, for approximately 2 seconds.
 - Pressure graphic is displayed.
- Switch off by pressing on/off again
 - Display off

Select language and pressure unit

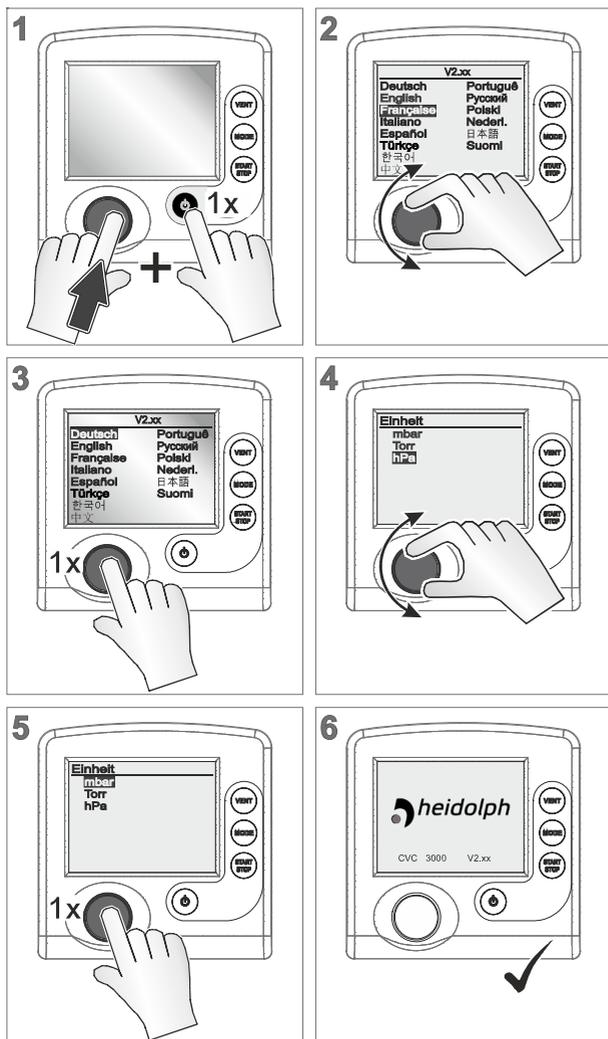Select operation mode

The controller is supplied with several operating modes. Only when the controlling process is stopped it is possible to select an operation mode.

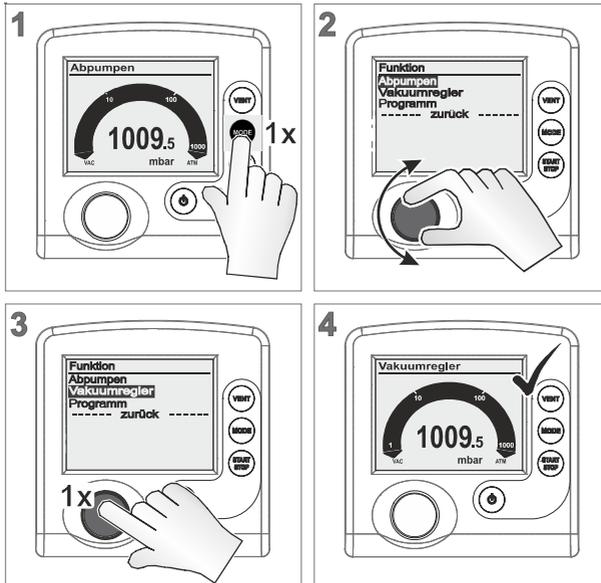Pump down

- Continuous pump down.
- Pump down with pressure and time presetting.
- VARIO®: pump down with adjustable motor speed (pumping speed) and continuous speed control..

Vac control

- Control to a preset vacuum value.

Program

- Load, edit and/or store program.
- Max. 10 programs with pressure and time presetting.


Select the mode suitable for vacuum apparatus and planned process.
Select any other operation mode in the same way as described above for Vac control.

After 20 seconds without action, the display will return automatically to pressure display..

Start controlling

Start vacuum controlling after selecting the required operation mode. The controller works in delivery status with the default settings of the factory setting.

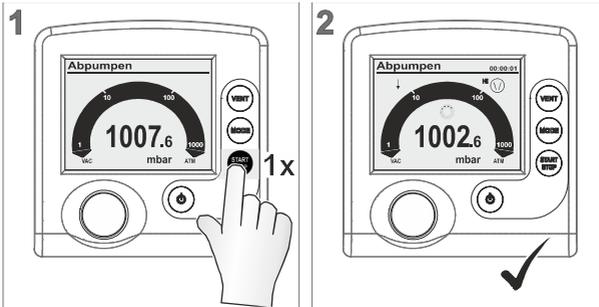- Controller starts.
- Icons are displayed.

Venting

DANGER

Danger of explosion when venting with air by forming of explosive mixtures. Depending on the process venting can cause formation explosive mixtures.

- Never vent processes with air which can form explosive mixtures.
- If necessary vent with inert gas (max. 1.2 bar absolute).

Certain processes may cause overpressure.

The [VENT] button is used to vent the system. A short click on this button will momentarily vent the system as the process continues. Holding the [VENT] button for longer than 2 seconds will cause the system to be vented to atmospheric pressure and the pump will stop running; max. 1060 mbar (795 Torr). Continuous venting stops when pressing [VENT] button again.


- Momentarily venting: Venting impulse, venting valve respectively air admittance valve opens momentarily → short-term pressure increase.
- Continuous venting: Icon for venting valve is flashing, venting valve opens → continuous pressure increase until atmospheric pressure → venting valve closes. Controller stops..

Change operation mode

During running operation the operation mode can be switched between Pump down and Vac control by pressing Mode button.


Pump down → Vac control

Semiautomatic distillation. Recommended for applications for which the process vacuum is still to be determined. Firstly the vacuum pump is pumping down rapidly in mode Pump down. As soon as the required process vacuum has been reached, e. g., boiling vacuum, this vacuum can be maintained by switching to Vac control. The actual pressure is adopted as the required set vacuum.

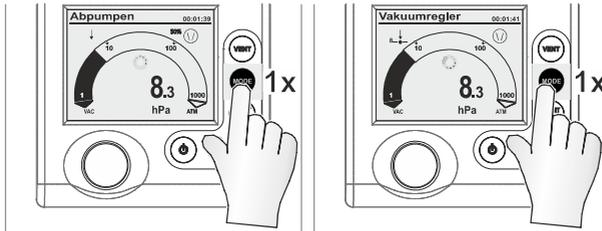Vac control → Pump down

With a connected VARIO® pump a controller working in Auto mode will detect and track the boiling point automatically. The vacuum will be adapted continuously to the process. If a particular process vacuum is required, the mode can be switched back again to Vac control.

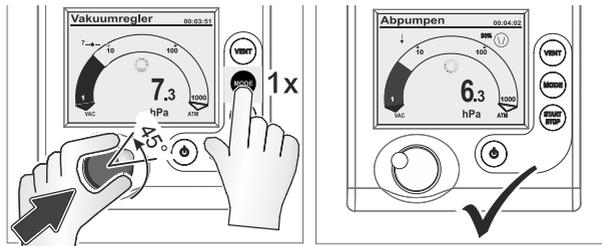
Switching the mode during operation via Mode button works only temporarily. After stopping the controller switches back to its primary mode.

Display graphic (curve)


In addition to the bar graphic of the pressure display the display can be switched to a diagram named Graphic which shows a pressure vs. time curve.

That Graphic curve will only be displayed while operation is running. With each start the recording restarts.

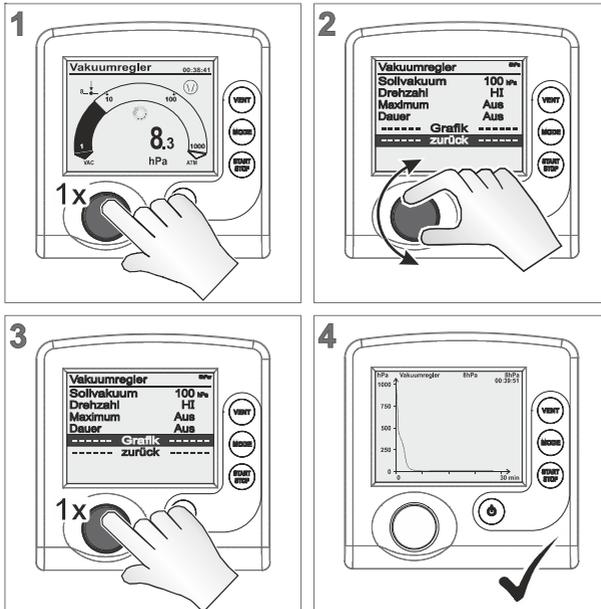To call up Graphic with pressure history for other operation modes, do like described above.


Example

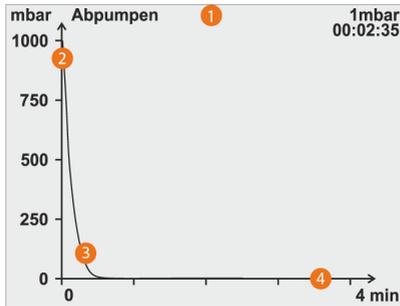Legend

Header	
1	Active operation mode Vacuum set point; for Vac control or Auto mode Actual vacuum value (actual pressure) Elapsed process time
Pressure axis	
2	Unit according to pre-settings (mbar, Torr, hPa)
Pressure graph	
3	Pressure/time progress
Time axis	
4	Continuous, automatically scaling time (minute, hour)

Quit display graphic

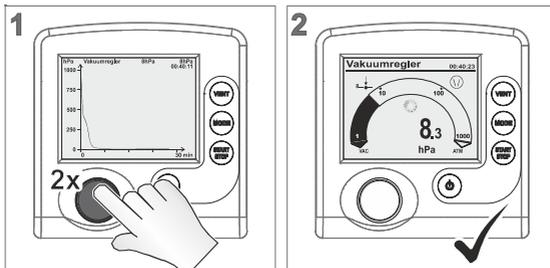
Set vacuum

In mode Vac control the set vacuum can be adapted directly during running operation.


Adapt set vacuum – quick tuning

- Press selection knob and turn it clockwise to increase set vacuum (venting).
- Press selection knob and turn it anticlockwise to decrease set vacuum (vacuum pump on).
 - Controller controls to the new set vacuum which is displayed while releasing the selection knob:

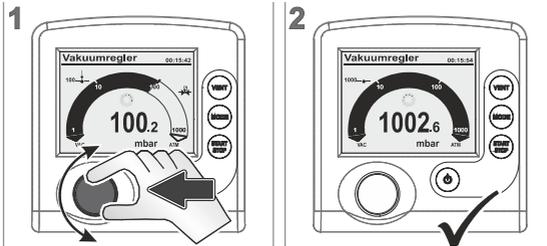
Adapt set vacuum – fine tuning

- 1 detent = 1 pressure value (mbar, Torr, hPa)
 - Controller regulates to new set vacuum:

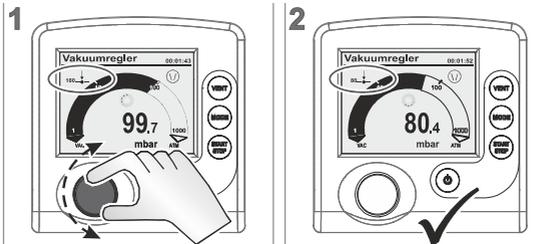Stopping the controller

Press [Start/Stop] to stop the controller.

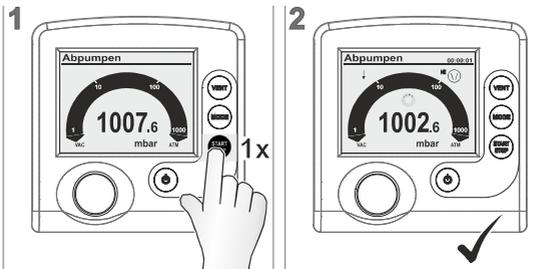
- Controller and vacuum control stops.
- Display icons switched off.

Operation menus

A selected operation mode can be adapted and optimized for the process through the corresponding operation menu. The settings in an operation menu include mainly: motor speed, set vacuum or time presettings. Settings in operation menus are retained also after switching on/off.

Menu Program is for storing up to 10 individual programs, e. g., to store control settings for frequently repeated processes.

Call up submenu of an individual operation mode

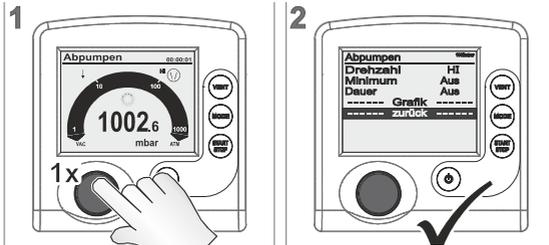To call up an operation menu, press the selection knob. Adaptions are possible during a running process as well as when control is stopped.

Use operation menu to optimize vacuum control for application requirements.

Pump down

Continuous pump down with pressure and time pre-settings..

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Parameter	Meaning
Minimum (mbar, Torr, hPa)	Vacuum set point; once reached, the controller switches off the vacuum pump or closes the in-line valve. Adjustment range: Off; 1-1060
Duration (min)	Pre-setting process runtime from Start on. Adjustment range: Aus; 1-1440

If Minimum and Duration are set to OFF , pump down has to be stopped by pressing START/ STOP button.

Application example – cabinet dryer

Set Minimum to a vacuum value below boiling pressure and the controller will switch the vacuum pump off, once the liquid has completely evaporated.

Vac control

Control to a set vacuum value..

Parameter	Meaning
Set vacuum (mbar)	Setting for lower vacuum level for 2-point control or precisely for VARIO® pump. Adjustment range: Turbo*; 1–1060
Hysteresis (mbar, Torr, hPa)	Only for VMS or in-line valve with vacuum pump: control range for 2-point control. Adjustment range: Auto; 1–300
Maximum (mbar, Torr, hPa)	Setting for upper vacuum level. Once reached, control switches off. Adjustment range: Off; 1–1060
Duration (min)	Presetting process runtime from Start on. Adjustment range: Off; 1–1440

Turbo mode: auto-adapting vacuum control for best ultimate vacuum. Best backing pressure for operation with a turbomolecular pump.

Application example – filtration


Set the set vacuum higher than the boiling pressure of the liquid and set Maximum value even a little bit higher. If the filter runs dry or if the filter is fractured, the pressure will increase and the control will be stopped automatically.

Hysteresis values Auto

Set vacuum (mbar)	5	10	50	80	100	200	500	700	900	1000
Hysteresis (mbar)	2	2	5	8	9	17	40	55	71	78

Program

Up to 10 individual programs including vacuum and time presettings can be loaded, edited, and stored.

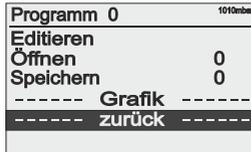Parameter	Meaning
Edit	Edit program with presettings for a process cycle or edit an existing program.
Open	Load the selected program.
Store	Store the program under the selected number. (memory capacity for up to 10 programs)
Hysteresis (mbar, Torr, hPa)	Only for VMS or in-line valve with vacuum pump: control range for 2-point control. Adjustment range: Auto; 1–300


For further descriptions of program functions, please contact an authorized dealer or our technical service, see „Contact information Heidolph international“ on page 102.

Hysteresis values Auto

Set vacuum (mbar)	5	10	50	80	100	200	500	700	900	1000
Hysteresis (mbar)	2	2	5	8	9	17	40	55	71	78

Configuration

In menu Configuration the controller parameters are set. This menu is also for adjusting the vacuum sensor and for loading Defaults settings.

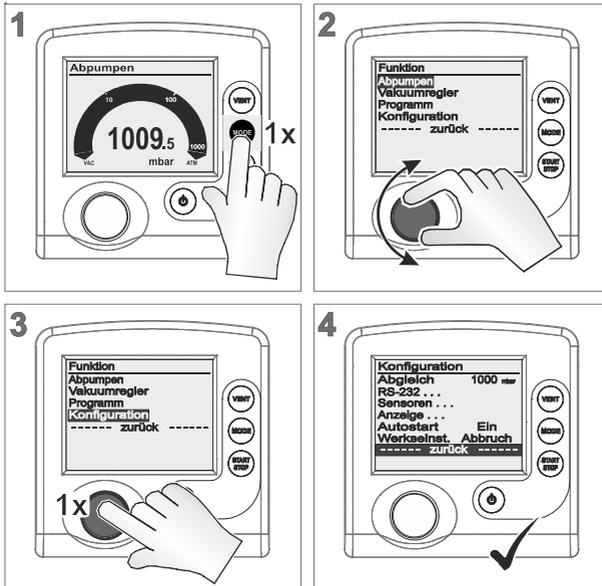After 20 seconds without action, the display will return automatically to pressure display.

Default settings

If Defaults Load is activated all controller parameter will be reset to delivery status. Stored programs will be deleted.

Options	Meaning
Load	Load default settings.
Cancel	Leave menu item without default setting.

Autostart

Ensure, if Autostart is activated, that no hazardous situations may occur due to the automatic start of the process.

Check whether the Autostart feature can be used safely with the intended application.

Options	Meaning
Off	After switching on power supply or after power failure the controller remains in Stop. Press Start/Stop button to start the controller.
On	Once power is applied, the controller starts automatically with the settings before power failure. The controller starts control directly without pressing Start/Stop button, if it previously was in running operation. Recommended, if power supply is switched on from a central point or if power is switched on by an external switch.

Display

Options	Meaning
Brightness	Adjust backlight brightness of the display. (0–100 %).
Contrast	Adjust display contrast. (0–100 %).
Sound	Off – Switch off keystroke sound and warning sound. On – Switch off keystroke sound and warning sound.
Unit	Preset pressure unit for user interface.
Language	User interface language

Sensors

In submenu Sensors all connected sensors are listed. The internal sensor is generally displayed as Sensor. External sensors are listed with sensor type name and address.

Sensoren	
Sensor	990.8 mbar
VSP 1	4.1E+2 mbar

Options	Meaning
Inverse	Sensor = currently selected for pressure display.
Sensor type	Selection for displaying pressure on basic display (max. 8 sensors are listed).

If required, contact an authorized dealer or our technical service, see „Contact information Heidolph international“ on page 102.
The display switches automatically to the previous menu when selecting a sensor with the selection knob.

RS-232

Submenu RS-232 is applied for interface configuration, parameter adjustments and commands.

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Parameter	Options	Meaning
Baud	<ul style="list-style-type: none"> ▪ 19200 ▪ 9600 ▪ 4800 ▪ 2400 	<p>Default setting for transmission speed.</p> <p>The baud rate of data transfer of transmitter and receiver must correspond.</p>
Parität	<ul style="list-style-type: none"> ▪ 8-N-1 ▪ 7-O-1 ▪ 7-E-1 	<p>Default setting for parity check, a method for error detection.</p>
Handshake	<ul style="list-style-type: none"> ▪ RTS-CTS ▪ Xon-Xoff ▪ none 	<p>Preset for continuous data transmission without loss – flow control.</p>
Remote	<ul style="list-style-type: none"> ▪ off 	<p>Control commands not enabled, only queries possible.</p>
	<ul style="list-style-type: none"> ▪ on 	<p>Connection for communication via RS 232 interface enabled.</p>

When selecting Remote On the controller itself is only operable via an external device. All keys of the control panel except button On/Off are locked.

VACUU-CONTROL® detects automatically, if Remote is activated or deactivated and retains that setting.

Function

If required, please contact our technical service, see „Contact information Heidolph international“ on page 102.

Program functions

If detailed descriptions about function menu with address assignment or about differential pressure measurement are required, please the technical service of Heidolph Instruments.


Troubleshooting

Error display

The major symbol for fault indication is the warning triangle. Additionally displayed icons and sounds refer to the cause of fault:

	Flashing: Warning where applicable with <ul style="list-style-type: none"> flashing component icon, warning sound (only when switched on) or flashing backlight.
	in combination with number = Vacubus address of the component which is defective.

Display icons

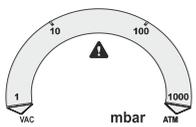Icon	Fault	Sound
	Limit pressure reached	1x
	Overpressure	1x
	Process time elapsed	1x
	Venting valve	2x
	In-line suction valve	3x
	Coolant valve	4x
	External sensor removed or defective internal sensor defective	5x 7x
	Vario pump	6x
	VACUULAN process pressure not reached within 99 hours	8x

	Digital I/O module: <ul style="list-style-type: none"> Fault indicator triggered Fault special configurations 	9x
	Level sensor triggered; flask full	10x
	Emission condenser Peltronic (too hot)	11x
	Analog I/O module	12x

A defective I/O module, which is configured as a remote module, does not trigger a warning alert. The control is stopped. Alert display by the flashing warning triangle.

Error messages

If a fault cannot be rectified with the described suggestions, please contact an authorized sales representative or our technical service (see „Contact information Heidolph international“ on page 102).

	possible cause	remedy
Sensitive process not controllable	Motor speed too high. Pumping speed too high.	<ul style="list-style-type: none"> Reduce motor speed.
Frequent error messages of connected components	Several controllers are connected. Several VACUU-BUS components of the same type are using the same address.	<ul style="list-style-type: none"> Use only one controller for one VACUU-BUS system. Only in Function menu Assign a new address number to VACUU-BUS@ component
VENT button does not work	Venting function deactivated. Setting of int. Air V in menu Function is switched Off or Auto.	<ul style="list-style-type: none"> Check why Venting is deactivated. Check if Venting by internal air admittance valve can be used without risk.
Internal air admittance valve cannot be triggered	External air admittance valve is connected, and/or External vacuum sensor is connected.	<ul style="list-style-type: none"> Venting safe? Enable the function in menu Function/ int. Air V; adjustment: Auto or On.
Internal air admittance valve does not switch	Air admittance valve soiled.	<ul style="list-style-type: none"> Clean venting valve (air admittance valve).
Function or menu item cannot be used	Function or menu item possibly only usable with short-cut (button combination).	<ul style="list-style-type: none"> Press the correct button combination

Vario pump icon flashes	<p>VARIO pump and VMS are both connected at the same time.</p> <p>VARIO pump defective.</p> <p>VMS defective or cable is not conneted.</p> <p>Cable break.</p>	<ul style="list-style-type: none"> ▪ Remove VMS from VARIO-pump and restart controller. ▪ Check VARIO pump for defective parts. ▪ Check VMS for defective parts. ▪ Check cable connections. ▪ Replace defective parts.
Air admittance valve icon flashes	<p>External air admittance valve removed.</p> <p>Plug disconnected.</p> <p>External air admittance valve defective.</p>	<ul style="list-style-type: none"> ▪ Check the connection. ▪ Check plug connection. ▪ Replace defective parts. ▪ Use internal air admittance valve. ▪ Reconfiguration without air admittance valve.
In-line suction valve icon flashes	<p>In-line suction valve removed.</p> <p>Plug disconnected.</p> <p>In-line suction valve defective.</p>	<ul style="list-style-type: none"> ▪ Check the connection. ▪ Check plug connection. ▪ Replace defective parts. ▪ Reconfiguration without In-line suction valve. ▪ Switch-off the controller; On/Off button. ▪ Remove In-line suction valve and switch on controller again.
Coolant valve icon flashes	<p>Coolant valve removed.</p> <p>Coolant valve defective.</p>	<ul style="list-style-type: none"> ▪ Check the connection. ▪ Replace defective parts. ▪ Reconfiguration without coolant valve.
Level sensor icon flashes	<p>Level sensor triggered (flask full).</p> <p>Level sensor removed.</p> <p>Level sensor triggered with empty flask.</p> <p>Cable break.</p> <p>Level sensor defective.</p>	<ul style="list-style-type: none"> ▪ Empty flask/catch pot. ▪ Check position of level sensor. ▪ Adjust level sensor or delete the sensor from controller (by loading default). ▪ Check plug connection. ▪ Replace defective parts.
Peltronic icon flashes	<p>Peltronic emission condenser too hot.</p> <p>Plug disconnected.</p>	<ul style="list-style-type: none"> ▪ Let the Peltronic emission condenser cool down. ▪ Check plug connection.
Title bar without text	<p>No controllable device connected (In-line suction valve, VMS, VARIO pump).</p>	<ul style="list-style-type: none"> ▪ Check device connections and cable. ▪ Replace defective parts. ▪ Connect a controllable device to the controller. ▪ Use the controller as measuring gauge.

Troubleshooting

No button reaction – only On/Off, PC icon displayed	Remote switched On. Controller only controllable via connected external end device (via RS232).	<ul style="list-style-type: none"> ▪ Switch-off Remote (switch off and on cotroller, press selection knob shortly while booting, select Configuration/RS232/Remote and adjust Off). ▪ Control controller via end device.
No reaction to button actuation	Controller defective	<ul style="list-style-type: none"> ▪ Contact us and ▪ return device for repair.
No display	Controller switched off. Power supply disconnected. Power supply not correctly connected. Mains voltage failure. Controller defective Cable break.	<ul style="list-style-type: none"> ▪ Switch on the controller; On/Off button. ▪ Check plug connection and wall power supply for correct connection. ▪ Replace defective parts. ▪ Contact service and ▪ return device for repair.
Blank display	Too many devices connected, e. g., valves. Short circuit of a connected device. Short circuit at RS232 interface. Controller defective	<ul style="list-style-type: none"> ▪ Power input of all connected devices may not exceed the maximum power consumption of the controller: <ul style="list-style-type: none"> - controller with wall power supply max. 30 W, - Controller + VARIO max. 25 W. ▪ Replace defective parts. ▪ Check RS232 plug connection. ▪ Contact Service and ▪ return device for repair.
Incorrect pressure display	Humidity inside the vacuum sensor. Vacuum sensor soiled. Vacuum sensor not adjusted. Vacuum sensor not correctly adjusted.	<ul style="list-style-type: none"> ▪ Identify and remove source of humidity. ▪ Dry the vacuum sensor, e. g., by pumping down. ▪ Clean the vacuum sensor. ▪ Readjust vacuum sensor.
Digital pressure gauge flashes	Pressure display flashing with O.O: - vacuum adjustment not correctly carried out. Pressure display flashing: - Overpressure! Pressure > 1060 mbar.	<ul style="list-style-type: none"> ▪ Readjust internal or external vacuum sensor. ▪ Risk of bursting. ▪ Discharge the system immediately by venting.
No digital pressure reading	External vacuum sensor defective. External vacuum sensor removed. Internal vacuum sensor defective.	<ul style="list-style-type: none"> ▪ Replace defective parts. ▪ Reconnect external vacuum sensor. ▪ Contact Service and ▪ return device for repair.
Sensors submenu is permanently displayed	Submenu Sensors... does not automatically switch back to previous display.	<ul style="list-style-type: none"> ▪ Select the required sensor by turning and pressing selection knob.

After loading defaults Language selection appears	Special factory settings have been loaded.	<ul style="list-style-type: none"> Set language and pressure unit. IMPORTANT! Check if the loaded default settings are suitable for your vacuum apparatus.
Error I/O module	<p>Plug disconnected.</p> <p>An error occurred in the system, the I/O module passed the error alert to the controller.</p>	<ul style="list-style-type: none"> Check plug connection. Remedy external fault.
VSP sensor displays wrong values	VSP sensor configured as VSK.	<ul style="list-style-type: none"> Use menu Function/Vacuubus to reconfigure the sensor as VSP.
Controller in operation, pressure display flashes	VSK sensors are measuring negative difference pressure.	<ul style="list-style-type: none"> Select an other vacuum sensor in menu Sensors.
Elapsed process time	<p>All program steps are completed.</p> <p>Program end reached.</p>	<ul style="list-style-type: none"> Acknowledge indication by pressing Start/Stop button.
Flashing clock icon	Elapsed process time	<ul style="list-style-type: none"> Acknowledge indication by pressing Start/Stop button.
Pump down stops, flashing arrow down icon	Pressure below preset minimum value.	<ul style="list-style-type: none"> Acknowledge indication by pressing Start/Stop button. If possible readjust presetting (Min.).
Vac control stops, flashing arrow up icon	Preset maximum value exceeded.	<ul style="list-style-type: none"> Acknowledge indication by pressing Start/Stop button. If possible readjust presetting (Max.).
Program -	Program not yet stored.	<ul style="list-style-type: none"> Store program under a free program number.

Error messages external components

Error messages for defective external components such as Inline suction valve, vacuum sensor, etc. cannot be reset. Replace defective accessories or send defective accessories for repair to your local supplier or to our service.

Reset

Auto reset

The following error indications will be reset automatically with remedy:

- Overpressure
- Process time elapsed
- Limit pressure reached
- Error air admittance valve
- Error Peltronic

Active reset

Several error indications need to be reset manually. Depending on the fault severity different actions are required.

- Press Start/Stop button to reset the following error indications:
 - In-line suction valve error
 - Coolant valve error
 - External vacuum sensor removed
 - I/O module activated Error indication
 - external error indicator has triggered via Digital I/O module; assigned as Error.
 - Level sensor triggered
- Load Defaults (standard factory setting) to reset the following error indications:
 - Missing set value presetting or VACUU-BUS plug disconnected via Digital I/O module; assigned as Remote.
 - Level sensor removed and/or VACUU-BUS plug disconnected.

Note the settings of stored programs, before loading Default.

Technical data

General device data	
Model	Vac control Hei-VAC Control
Dimensions basic unit (w × h × d)	123 × 124 × 83 mm (5 × 5 × 3.5 in.)
Dimensions with foot (w × h × d)	144 × 124 × 115 mm (6 × 5 × 4.5 in.)
Weight basic unit	440 g (0.97 lb)
Weight with foot	570 g (1.3 lb)
Display	LC, with brightness control
Pressure sensor	Integrated (ceramic diaphragm, alumina, capacitive, gas independent, absolute pressure)
Degree of protection (front)	IP42 (IP42)
Electrical data	
Rated voltage	24 VDC (±10 %)
Power input (max.)	3,4 W
Max. admissible current total for connected valves	4 A
Interface	RS232, SUB-D, 9pol
Permissible ambient conditions	
Operating temperature	10 °C – 40 °C (50 °F – 104 °F)
Transport/storage temperature	-10 °C – 60 °C (14 °F – 140 °F)
Rel. humidity	30 – 85 %, w/o condensation
Installation altitude	up to 3,000 m (9,840 ft) asl
Vacuum data	
Measuring rang, absolute	1080 – 0,1 mbar (810 – 0.1 Torr.)
Control range	1060 – 0,1 mbar (795 – 0.1 Torr.)
Resolution	0,1 mbar (0.1 Torr.)
Vacuum sensor VSK 3000	
Max. admissible pressure, abs.	1,5 bar (1125 Torr.)
Max. admissible media temperature (gas):	Temporarily: 80 °C (176 °F) Continuous: 40 °C (104 °F)
Measurement uncertainty	<±1 mbar (<±0.75 Torr.)
Temperature coefficient	<±0,07 mbar/K (<±0.05 Torr./K)
Venting	
Max. admissible pressure, abs.	1,2 bar (900 Torr.)
Connection	Hose nozzle for flexible tube inner d = 4 –5 mm

Wetted parts


Component	Material
Vacuum connection, hose nozzle	PP
Sensor	Aluminium oxide ceramic
Sensor housing	PPS/Glasfaser
Sensor seal	chemically resistand fluoroelastomer
Venting valve seal	FFKM

Interface commands

If required, contact an authorized dealer or our technical service, see „Contact information Heidolph international“ on page 102.

Pin assignment RS232 interface

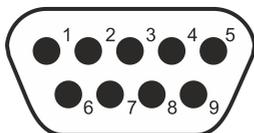Pin	Denomination	Function
1	DCD	–
2	RxD	Received data
3	TxD	Transmitted data
4	DTR	+10 V
5	GND	Ground
6	DSR	–
7	RTS	Transmission request
8	CTS	Ready to send
9	RI	+5 V (Bluetooth, remote control)

Scope of delivery

Item	Quantity	Product no.
Controller Hei-VAC Control	1	11-001-591-28
Power supply unit 30W 24V; including interchangeable mains plugs	1	11-300-008-12
Vacuum valve preassembled	1	569-00080-00
Plastic foot CVC 3000	1	11-300-008-11
User manual	1	01-005-005-87
Warranty registration	1	01-006-002-78

Accessories

Detailed information on the available accessories for your device variant can be found on our website at www.heidolph-instruments.com.

In case of need, contact an authorized dealer or our technical service, see „Contact information Heidolph international“ on page 102.

Device service

WARNING

Live components are installed inside the device. When opening the device, there is a risk of touching live components.

- Disconnect the device from the power supply before carrying out any cleaning, maintenance or repair work.
- Penetrating liquid poses the risk of an electric shock! When cleaning, avoid the penetration of liquids.

General cleaning instructions

Wipe all surfaces of the device with a damp cloth if necessary. Persistent contamination can be removed with mild soapy water.

CAUTION

In the event of contact with aggressive media, make sure that all surfaces are cleaned immediately with suitable aids to prevent damage to the device surface.

Improper cleaning can damage the surfaces of the device. Furthermore, penetrating liquid can damage the electronic components inside the device.

- Clean the device's surfaces with a soft, lint-free and only slightly moistened cloth.
- Never use any aggressive or abrasive cleaning agents or aids.

Clean venting valve

- Apply slight overpressure of dry air or inert gas to the vacuum port.
- Press the VENT key several times until gas escapes through the venting port.
- Repeat this procedure until you hear the clicking of the valve and a gas stream is noticeable at the venting port.

Clean internal sensor

- Fill a small amount of solvent via the vacuum port (1) in the controller, e. g., cleaning solvent.
- Let the solvent react for a few minutes.
- Drain the solvent.
- Dissolved substances or discolorations in the solvent are possible.
- Repeat this procedure until no more pollutants are in the solvent.
- Let the controller dry.
- Readjust the internal (vacuum) sensor.

Sensor adjustment

For readjustment the reference pressures need to be known with certainty. In the pressure range 20 – 700 mbar (15 – 525 Torr) no adjustment is possible.

Check the accuracy of the pressure sensor in case of irregularities in the pressure display. Adjust the sensor in two steps: at atmospheric pressure and under vacuum.

Do not adjust at atmospheric pressure, if the pressure at the location of the device is not exactly known (pay attention to height above sea level).

Any kind of pollution of the vacuum system, e. g., oil, substances, or humidity could falsify the adjustment. Clean polluted sensors before readjustment.

Adjustment at atmospheric pressure

An adjustment at atmospheric pressure is only possible if the pressure is higher than > 700 mbar (> 525 Torr).

- Vent the measurement connection of the controller or in case the connected external vacuum sensor VSK 3000.
- Make sure that the vacuum sensor (internal or external) is really at atmospheric pressure.
- Determine the exact atmospheric pressure of your location, e. g, by barometer, inquiry at the meteorological office or the airport.
- Call up menu Configuration.
- Turn the selection knob and place the bar marking on Adjustment.
- Press the selection knob.
 - Marking jumps to numeric value.
- Adjust the exactly determined local atmospheric pressure by turning the selection knob.
- Press the selection knob.
 - Sensor adjusted to atmospheric pressure.

Adjustment under vacuum

An adjustment under vacuum is only possible if the pressure is lower than < 20 mbar (< 15 Torr) absolute.

- Evacuate the measurement connection of the controller or in case the connected external vacuum sensor VSK 3000 to a pressure < 0,1 mbar.

Adjustment under vacuum with an actual pressure higher than 0,1 mbar (0.1 Torr) reduces the accuracy of the measurement. If the pressure is significantly higher than > 0,1 mbar (> 0.1 Torr) the adjustment to a reference pressure is recommended.

- Call up menu Configuration.
- Turn the selection knob and place the bar marking on Adjustment.
- Press the selection knob.
- Marking jumps to numeric value.
- Adjust the pressure value to 0 by turning the selection knob.
- Press the selection knob.
 - Sensor adjusted under vacuum.

The readjustment of a VSP 3000 can only be carried out in warmed-up state.

Adjustment is not possible during the warm-up time.

Use a high vacuum pump for the adjustment of a VSP sensor.

After connection to power supply and after the pressure has reached < 10-3 mbar, wait 20 minutes before adjusting the VSP sensor.

Carry out the adjustment in the same order as described above for VSK.

Adjustment at a reference pressure

Instead of adjustment under vacuum to a pressure < 0,1 mbar (< 0.1 Torr), adjustment to a precisely known reference pressure within the range of 0 – 20 mbar (0 – 15 Torr) is possible.

- Evacuate the measurement connection of the controller or in case the connected external vacuum sensor VSK 3000 to a pressure in the range of 0 – 20 mbar (0 – 15 Torr).
- Call up menu Configuration.
- Turn the selection knob and place the bar marking on Adjustment.
- Press the selection knob.
- Marking jumps to numeric value.
- Adjust the pressure value to the actual reference pressure by turning the selection knob.
- Press the selection knob.
- Sensor adjusted to reference pressure.

The measurement uncertainty of the reference pressure will directly affect the measurement uncertainty of the controller.

If the nominal ultimate vacuum of a diaphragm pump is used as reference vacuum, the accuracy of the controller might be doubtful. The diaphragm pump may not achieve the specified vacuum (due to condensate, poor condition, failure of valves or diaphragm, leaks).

Repairs

Repairs to the device may only be carried out by authorized experts! Unauthorized repairs during the warranty period will result in the loss of the warranty claim. The owner is solely liable for damage caused by unauthorized repairs.

If repairs are required, contact an authorized dealer or our technical service, see „Contact information Heidolph international“ on page 102. Include the completed certificate of decontamination with every device return, see „Certificate of decontamination“ on page 103.

Maintenance

There are no user-serviceable components in the unit housing. If necessary (in the event of abnormal operating behavior such as excessive noise or heat generation, for example), contact our technical service, see „Contact information Heidolph international“ on page 102.

Disposal

- When disposing of the device, observe the provisions of the WEEE Directive 2012/19/EU and its transposition into national law in the country of use.
- When disposing of portable batteries, observe the provisions of the European Battery Directive 2013/56/EU and their implementation in national law in the country of use.
- Check the device and all components for residues of substances that are hazardous to health, the environment and biohazardous before disposal.
- Properly remove and dispose of residues of substances that are hazardous to health, the environment and are biohazardous!

Contact information Heidolph international

Heidolph Instruments North America

Phone: 1-866-650-9604
 E-mail: service@heidolph.com
www.heidolphNA.com

Heidolph Instruments United Kingdom

Phone: 01799 - 5133-20
 E-mail: service@radleys.co.uk
www.heidolph-instruments.co.uk

Local distributors

To find your local distributor please visit www.heidolph.com

Warranty Statement

Heidolph Instruments provides a three-year warranty against material and manufacturing defects.

Glass and wear parts, transportation damage, and damage resulting from improper handling or non-intended use of the product are excluded from the warranty.

The warranty period for registered products begins on the date of purchase. Register the product with the enclosed warranty card or on our homepage www.heidolph.com.

For non-registered products, the warranty period begins with the date of the serial production (to be determined by the serial number).

In the event of material or manufacturing defects, the product will either be repaired or replaced free of charge within the warranty period.

Certificate of decontamination

Enclose the Certificate of decontamination, duly completed, with your device return.
Submissions without a Certificate of decontamination cannot be processed!

Certificate of decontamination

IN CASE OF RETURNS

Please fill in the required fields.

Note: The sender must package the goods properly and appropriately for transport.

Heidolph Instruments GmbH & Co. KG
Walpersdorfer Straße 12
91126 Schwabach, Germany
Phone: +49 (0) 9122 9920-380
Fax: +49 (0) 9122 9920-19
Email: service@heidolph.de

SENDER

Name _____ First name _____
Company/Institution _____ Department _____
Workgroup _____
Address _____
ZC/City _____
Country _____ Phone _____
Email _____

DEVICE DETAILS

Article number _____ Serial no. _____
Ticket number _____
Reason for sending in _____

Has the device been cleaned, decontaminated/disinfected? Yes No (Please mark as applicable)
If yes, which measures were carried out?

Does this device pose a risk to people and/or the environment due to the processing of substances that are hazardous to health, the environment and/or are biohazardous? Yes No (Please mark as applicable)
If yes, with which substances did the device come into contact?

LEGALLY BINDING DECLARATION

The principal/consignor is aware that they are liable to the agent/consignee for losses or damage incurred due to incomplete and incorrect information.

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hereby declare, that the product (from serial number on 200150525) designated below is in compliance with the basic requirements of all applicable EU-directives stated below with regard to design, type of model sold and manufactured by us. This certificate will be invalid if the product is modified without the prior written consent and agreement of the manufacturer.

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UKCA Declaration of Conformity

**UK
CA**

Declaration of Conformity In accordance with UK Government guidance

This declaration is issued under the sole responsibility of the manufacturer,
Heidolph Instruments GmbH & Co. KG
Walpersdorfer Straße 12
91126 Schwabach / Germany

Product: Laboratory vacuum controller

Model:

Hei-VAC Control 591-00360-00

Description:

Vacuum controller for rotary evaporators

The object of the declaration described above is in conformity with the relevant
UK Statutory Instruments (and their amendments):

2016 No. 1101

2008 2016 No. 1091

2012 No. 3032

Electrical Equipment (Safety) Regulations 2016

The Electromagnetic Compatibility Regulations 2016

*The Restriction of the Use of Certain Hazardous Substances in
Electrical and Electronic Equipment Regulations 2012*

and complies with the following technical standards :

EN ISO 12100:2010, EN 61326-1:2013, EN 61010-1:2010, EN IEC 63000:2018

UK Authorised Representative (for authorities only):

ProductIP (UK) Ltd.
8. Northumberland Av.
London WC2N 5BY

Signed for and on behalf of Heidolph Instruments GmbH & Co. KG
Walpersdorfer Straße 12, 91126 Schwabach / Germany

Schwabach, 01.09.2021

Wolfgang Jaenicke
Managing Director

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With this declaration, we confirm (according to current knowledge) that all sold laboratory devices by Heidolph Instruments GmbH + Co. KG fulfill the requirements of the EU directive 2011 / 65 / EU (RoHS) and the amended of directive 2015 / 863.

All devices are compatible with the requirement of the RoHS for the following materials:

Max. 0,01% of the weight in homogeneous material for cadmium and max. 0,1 % of the weight in homogeneous material for lead, mercury, hexavalent chromium, polybrominated biphenyl, polybrominated diphenyl ether, Di (2-ethylhexyl l) phthalate, butyl benzyl phthalate, dibutyl phthalate, diisobutyl phthalate.

In the case of individual assemblies, maximum concentrations maybe exceeded within the permissible exceptions of the Directive.

Schwabach, 22.07.2019

Stefan Peters
Vice President Marketing, Innovation & Technology

Marcell Sarré
Vice President Quality Management & Technical Service

China RoHS Certification

China RoHS DECLARATION OF CONFORMITY

Heidolph Instruments GmbH & Co.KG has made reasonable efforts to ensure that hazardous materials and substances may not be used in its products.

In order to determine the concentration of hazardous substances in all homogeneous materials of the subassemblies, a "Product Conformity Assessment" (PCA) procedure was performed. As defined in GB/T 26572 the "Maximum Concentration Value" limits (MCV) apply to these restricted substances:

- Lead (Pb): 0.1%
- Mercury (Hg): 0.1%
- Cadmium (Cd): 0.01%
- Hexavalent chromium (Cr(VI)): 0.1%
- Polybrominated biphenyls (PBB): 0.1%
- Polybrominated diphenyl ether (PBDE): 0.1%

Environmental Friendly Use Period (EFUP)

EFUP defines the period in years during which the hazardous substances contained in electrical and electronic products will not leak or mutate under normal operating conditions. During normal use by the user such electrical and electronic products will not result in serious environmental pollution, cause serious bodily injury or damage to the user's assets.

The Environmental Friendly Use Period for Heidolph Instruments GmbH & Co.KG products is 25 years.

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This table is created according to SJ/T 11364-2014

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Note: Table applies to all products. Some of the components or parts listed above may not be part of the enclosed product.

- O: 表示该有毒有害物质在该部件所有均质材料中的含量均在GB/T 26572规定的限量要求以下。
- O: Indicates that the above mentioned hazardous substance contained in all homogeneous materials of the part is below the required limit as defined in GB/T 26572.

- X: 表示该有毒有害物质至少在该部件某一均质材料中的含量超出GB/T 26572规定的限量要求。
- X: Indicates that the above mentioned hazardous substance contained in at least one of the homogeneous materials of this part is above the required limit as defined in GB/T 26572.

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Apart from the disclosures in the above table, the subassemblies are not intentionally manufactured or formulated with lead (Pb), mercury (Hg), cadmium (Cd), hexavalent chromium (CrVI), polybrominated biphenyls (PBB), and polybrominated diphenyl ethers (PBDE).

Products manufactured by Heidolph Instruments GmbH & Co.KG may enter into further devices or can be used together with other appliances.

With these products and appliances in particular, Heidolph Instruments GmbH & Co.KG will not take responsibility for the EFUP of those products and appliances.

Schwabach, 06.08.2021

Wolfgang Jaenicke
Chief Executive Officer CEO

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Subject to change without notice. The printed version of this document is not regularly updated. The latest issue of this document can be found by visiting our homepage.

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MATERIAL CONTENT DECLARATION FOR Heidolph Instruments GmbH & Co. KG PRODUCTS
© Heidolph Instruments GmbH & Co. KG